

Peripheral Balloon Expandable Stent System (QBX) 7 F

TECHNICAL SPECIFICATIONS

Description	Expandable Stent System Semi-Compliant OTW
Material	Cobalt Chromium alloy
Lengths	17 26 36 57 mm
Recommended Guidewire	0.035" (0.89 mm)
Recommended Introducer Sheath	7 F
Nominal Pressure	7 bar
Rated Burst Pressure (RBP)	10 bar
Working Catheter Length	80 cm 120 cm
X-ray Marker	One marker at proximal and distal end of balloon
Entry Tip Profile	min 0.95 mm (0.037")
Crossover Profile	max 2.45 mm (0.096")
Foreshortening	As low as 1%
Recoil	As low as 1%
Balloon Material	Pebax
Hydrophilic Coating	400 mm distal

ORDER INFORMATION

Usable catheter length **80 cm**

Diameter (mm)	Stent Length (mm) Balloon Length (mm)			
	17 20	26 30	36 40	57 60
8	08 QBX 08017	08 QBX 08026	08 QBX 08036	08 QBX 08057
9	08 QBX 09017	08 QBX 09026	08 QBX 09036	08 QBX 09057
10	08 QBX 10017	08 QBX 10026	08 QBX 10036	08 QBX 10057

Usable catheter length **120 cm**

Diameter (mm)	Stent Length (mm) Balloon Length (mm)			
	17 20	26 30	36 40	57 60
8	12 QBX 08017	12 QBX 08026	12 QBX 08036	12 QBX 08057
9	12 QBX 09017	12 QBX 09026	12 QBX 09036	12 QBX 09057
10	12 QBX 10017	12 QBX 10026	12 QBX 10036	12 QBX 10057

CE 0124




QualiMed Innovative Medizinprodukte GmbH
 Boschstraße 16, 21423 Winsen, Germany
 Phone +49 4171 6578 0
 Fax +49 4171 6578 11
 Web www.qualimed.de

Peripheral Vascular



Peripheral Balloon Expandable Stent System (QBX) 7 F

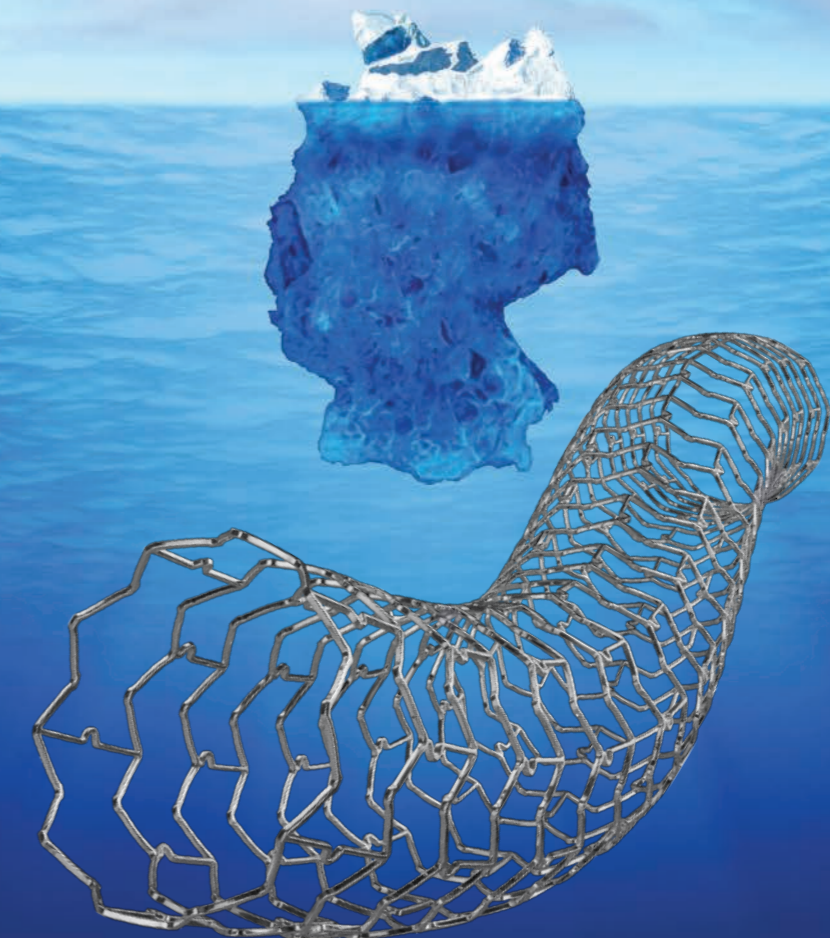

A flexible balloon-expandable stent for easy deployment and precise placement

Regulatory Status
CE Mark Approved



Peripheral Balloon Expandable Stent System (QBX) 7 F

The QBX 7 F is a flexible, balloon expandable stent, made of cobalt chromium alloy.

It is designed for peripheral vessel diameters from 8.0 mm to 10.0 mm. The QBX 7 F is mounted on an 0.035 over-the-wire delivery system and is available in a full range of diameters and lengths.

KEY FEATURES

- > **Low crossing profile** for all diameters
- > **Combination low shaft and tip entry profile**
- > **Soft, tapered balloon tip** for less trauma
- > **Flexible stent design**
- > **High radial strength**

BENEFITS

- > **Simple to use**
- > **Easy deployment and precise placement**
- > **Even and accurate stent expansion**
- > **Good accessibility and high pushability** across the peripheral vessels

COMPLIANCE TABLE

Pressure (bar)	Stent Diameter (mm)		
	8	9	10
5	7.76	8.8	9.64
6	7.88	8.9	9.82
7 NP	8.00	9.00	10.00
8	8.12	9.10	10.18
9	8.24	9.20	10.36
10 RBP	8.36	9.30	10.54
Nominal pressure	7	7	7
RBP	10	10	10

The Peripheral Balloon Expandable Stent System provides a means of safely advancing the stent through the peripheral vessels to the desired location.

